

Heel Snuggler™

Warmth Where You Need It

PEOPLE. PRODUCTS. PROGRAMS.®



Children's
Medical
Ventures

RESPIRONICS®

Children's Medical Venture's Heel Snuggler™ infant heel warmer wraps and rewraps securely around a baby's heel to offer warmth and comfort prior to a heel stick procedure. The unique butterfly design covers the sole and heel of the baby's foot, while the internal disk activates easily and provides warmth quickly and accurately. **Warmth where you need it.**

The Heel Snuggler™ is the only infant heel warmer designed to truly conform to the infant's heel and hold it firmly in place. The unique design is shaped like a butterfly and has two soft straps for a comfortable and secure fit. This new design provides consistent warming of the entire sole and heel of the foot prior to a heel stick procedure, as recommended by the NANN skin care guidelines, and removal is easy with perforated tabs.

Effective Warming:

- Unique design ensures complete coverage of heel stick area
- May be used any time increased capillary action is necessary
- Approximate temperature 105°F (41°C)



Specifications:

- Size: approximately 5.5" x 3.5" (14 cm x 9 cm)
- Ingredients: non-toxic sodium acetate (food grade), water and thickeners
- Temperature of the Heel Snuggler™ rises from 75°F (24°C) to 105°F (41°C) upon activation and remains at that temperature for three to four minutes

Safe and Easy to Use:

- Dual tapes reliably hold Heel Snuggler™ in place and can be retaped for perfect fit
- Straps "tear away" for easy removal
- Patented button activator eliminates risk of bag puncture
- Non-toxic ingredients
- Latex free

ordering information


item no.	description	quantity
99047	Heel Snuggler™	100/case 4 boxes of 25



Customer Service: 1-800-345-6443 or 724-387-4000
 Respironics Europe: +33-1-55-60-19-80
 Respironics Asia Pacific: +81-3-5800-0724

Heel Snuggler is a trademark of Respironics, Inc. and its affiliates. People Products. Programs. is a registered trademark of Respironics, Inc. and its affiliates.
 © 2004 Respironics, Inc. and its affiliates. All right reserved.

1023321 KB 9/27/04

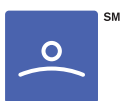

SLEEP
MANAGEMENT



CHRONIC RESPIRATORY
MANAGEMENT



TOTAL VENTILATION
SOLUTIONS



RESPIRATORY
DRUG DELIVERY



NEONATAL &
INFANT CARE