

SmartTrace™ Neonatal/ Pediatric Hospital Electrodes

The Smart, Skin-Friendly Choice



Children's
Medical
Ventures

RESPIRONICS®

Children's Medical Ventures' SmartTrace™ Neonatal/Pediatric Electrodes are designed with a baby's delicate skin, and condition, in mind. The electrode's unique construction adheres - and painlessly re-adheres - to the baby to reliably transmit critical information. This makes **SmartTrace™ Electrodes the smart, skin-friendly choice.**

SmartTrace™ Hospital Electrodes


SmartTrace™ Neonatal/Pediatric Hospital Electrodes are the smart, skin-friendly choice when you are looking for an electrode that stays sticky longer.

The electrodes are formulated with two unique patented gel layers. The first layer is extremely gentle on infant's skin and allows the electrode to stick and re-stick without harming the baby's delicate skin. The second layer is formulated for optimal adhesion to the electrode. These two layers are locked together with a unique mesh layer to prevent slippage/delamination.

- Patented dual layer hydrogel stays sticky longer
- Breathable moisture barrier performs well in high humidity
- Last up to 5-7 days with repeated application
- Double soldered wires reduce wire pullout
- Copper lead wire promotes higher quality tracings
- Applies and reapplies without rewetting
- Resists moisture and 'dry out'

Specifications:

- 0.06 in/2.05mm socket
- Pre-wired (black/white/green)
- Packaged 3 electrodes/set; 10 sets/case
- Two year shelf life

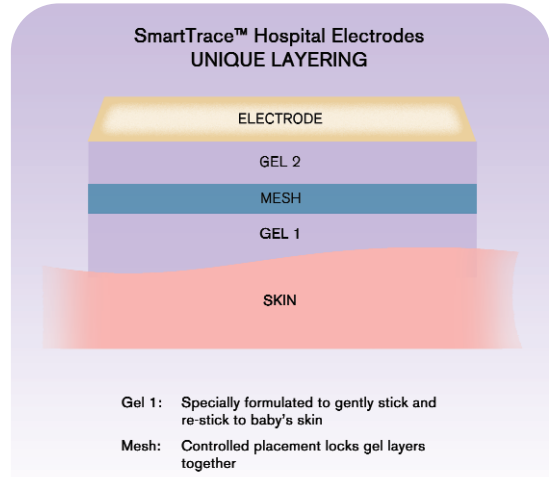

ordering information

item no.	description	quantity
1015663	SmartTrace™ Hospital Electrodes	10 sets/case
1015662	SmartTrace™ Homecare Electrodes	20 pairs/case



Customer Service: 1-800-345-6443 or 724-387-4000
 Respiroics Europe: +33-1-55-60-19-80
 Respiroics Asia Pacific: +81-3-5280-9611



SmartMonitor® 2 PS
Product Number 1014557

Use SmartTrace™ Hospital Electrodes with the SmartMonitor® 2: Professional Series or with any standard hospital monitor.

SmartTrace is a trademark of Respiroics, Inc. and its affiliates. SmartMonitor and People. Products. Programs. are registered trademarks of Respiroics, Inc. and its affiliates.
 ©2005 Respiroics, Inc. and its affiliates. All rights reserved.

1025739 (mci 4100139) 1/11/05 KB G

