

Wallaby® 3 Phototherapy System

Ideal for hospital and home use.



Children's
Medical
Ventures

RESPIRONICS®

Children's Medical Ventures' Wallaby® 3 Phototherapy System allows a mother to hold, even nurse, her baby without interrupting therapy. In fact, at just 7.2 pounds, it's light enough that it can be easily carried around with the baby. That's a much more natural way to deliver the precise level of therapy - and it gives new mothers more time to bond with their babies.

Wallaby® 3 Phototherapy System

The Wallaby® 3 Phototherapy System is a lightweight and compact alternative for delivering effective therapy to jaundiced babies. Every aspect of this device has been designed to make it the most complete - and most portable - phototherapy system available.

Efficient

Fiber Optic Technology:

- Softer, brighter pattern of therapeutic light - increased irradiance improved therapy delivery

Safe

Hour Meter:

- Automatically records treatment time.

Enhanced Safety Features:

- Locking device prevents the panel from disengaging from the illuminator if the cable is pulled.
- Light shuts off automatically if fiberoptic cable is disconnected.
- No infrared or ultraviolet light; not harmful to baby's eyes.

Flexible

Two Irradiance Levels:

- Increased light intensity to reduce higher bilirubin levels faster.

Two-Lamp Carousel with Automatic Bulb Exchange:

- Automatically exchanges the expired bulb, eliminating disruption of therapy and after-hours calls. A visual indicator alerts the caregiver that the second bulb is in use.

Two Sizes of Fiberoptic Panels:

- Choose the optimal modality for the infant: wrapped around the baby's torso or placed underneath the baby's body.



This neonatal phototherapy panel and the larger 14 inch wrap-around panel give you options to better tailor therapy to an individual baby's needs.



specifications

Illuminator Size	10.5" x 10.5" x 4" (26 cm x 26 cm x 10 cm)	
Weight	7.2 lbs (3.3 kg)	
Fiberoptic Panels	Neonatal panel: 4" x 5" (10 cm x 13 cm) Wrap-around panel: 3" x 14" (7.5 cm x 35.5 cm)	
Irradiance Levels	LEVEL 1	LEVEL 2
Neonatal panel	25-35 $\mu\text{W}/\text{cm}^2/\text{nm}$	50-60 $\mu\text{W}/\text{cm}^2/\text{nm}$
Wrap-around panel	10-20 $\mu\text{W}/\text{cm}^2/\text{nm}$	15-25 $\mu\text{W}/\text{cm}^2/\text{nm}$

ordering information


item no.	description	quantity
H-3600	Wallaby® 3 Phototherapy System	1/each
CC302	Wallaby® 3 carrying case	1/each
EG-2000	Wrap-around panel	1/each
EG-2000N	Neonatal panel	1/each
JD-101	Joey Dosimeter	1/each
MD-510	Disposable covers for wrap-around panel	50/pkg
MD-610	Disposable covers and vest for neonatal panel	100/pkg
PGM302	Parent's guide	1/each
HPK-301	Patient kit	1/each
H100-3600-20	Replacement phototherapy light bulb	1/each
161-5000-00	Replacement power cord	1/each



Customer Service: 1-800-345-6443 or 724-387-4000
 Respiroics Europe: +33-1-47-52-30-30
 Respiroics Asia Pacific: +81-3-5280-9611

Respiroics, the Children's Medical Ventures heart design logo and Wallaby are registered trademarks of Respiroics, Inc. and its affiliates.
 © 2005 Respiroics, Inc. and its affiliates. All rights reserved.

Gyer KB 10/10/05 MCI 4100375 PN 1029932

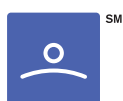

SLEEP
MANAGEMENT



CHRONIC RESPIRATORY
MANAGEMENT



TOTAL VENTILATION
SOLUTIONS



RESPIRATORY
DRUG DELIVERY



NEONATAL &
INFANT CARE