

# Enteral Only Extension Sets

When Safety Comes First



**Enteral Only Extension Set  
with Y-Site and Cap  
95017-C**

**Enteral Only Extension Set  
with Oral Adapter  
1017080**



Children's Medical Ventures' Enteral Only Extension Sets add a measure of safety to routine procedures. Our sets help reduce the possibility of connection errors that can occur when an infant is surrounded by lines, tubing and support equipment. The patented safety tip fits only into feeding tubes, not IV lines, helping to eliminate incorrect connections. The sets are also available with an oral adapter that accepts an oral dose syringe – another added safety feature and timesaving measure.

# Enteral Only Extension Sets

The Enteral Only Extension Sets are designed to help minimize the potential for inadvertent delivery of enteral feedings through the intravenous route. The unique safety tip connects easily to most naso-gastric (NG) tubes and enteral lines but is NOT compatible with IV tubing or stopcocks. Our Enteral Only Extension Sets are also available with an oral adapter that easily accepts an oral dose syringe directly from the pharmacy, helping to make administration of medications fast, easy and safe.







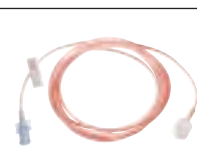
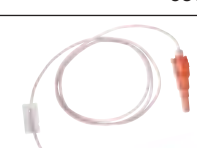

## Safety Features

- Safety tip connects only to feeding tubes and enteral lines
- Orange striping provides easy and quick identification of enteral feeding lines
- "Enteral Only" tag and slide clamp provide additional safety assurance

## Specifications

- Microbore PVC tubing
- Individual, sterile packaging
- Latex and DEHP-free

## Five Specialty Sets Available

	ENTERAL ONLY EXTENSION SETS	ENTERAL ONLY EXTENSION SETS WITH ORAL ADAPTER
<b>60" extension set</b> Priming volume: approx. 1.8 ml	 95017-A	 1017080
<b>60" extension set with stopcock</b> Priming volume: approx. 2.2 ml	 95017-B	 1018516
<b>60" extension set with Y-site and cap</b> Priming volume: approx. 2.5 ml	 95017-C	 1018517
<b>60" extension set with male luer end</b> Priming volume: approx. 1.8 ml	 95017-D	Not Available in This Configuration
<b>36" extension set</b> Priming volume: approx. 1.3 ml	 95017-E	 syringe not included 1018514

## ordering information

item no.	description	quantity
<b>Enteral Only Extension Sets</b>		
95017-A	60" extension set	50/cs
95017-B	60" extension set with stopcock	50/cs
95017-C	60" extension set with Y-site and cap	50/cs
95017-D	60" extension set with male luer end	50/cs
95017-E	36" extension set	50/cs
<b>Enteral Only Extension Sets with Oral Adapter</b>		
1017080	60" extension set	50/cs
1018516	60" extension set with stopcock	50/cs
1018517	60" extension set with Y-site and cap	50/cs
1018514	36" extension set	50/cs

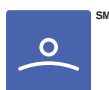


Customer Service: 1-800-345-6443 or 724-387-4000

Respironics Europe: +33-(0)1-55-60-19-80

Respironics Asia Pacific: +852-234-342-18

1021698 KB 5/24/04 G



SLEEP MANAGEMENT



CHRONIC RESPIRATORY MANAGEMENT



TOTAL VENTILATION SOLUTIONS



ASTHMA, ALLERGY & SINUSITIS MANAGEMENT



NEONATAL & INFANT CARE